


Doublecyte-N

Ago doppio lume per prelievo ovociti e lavaggio follicolare

Double lumen needle for oocyte aspiration and follicle flushing

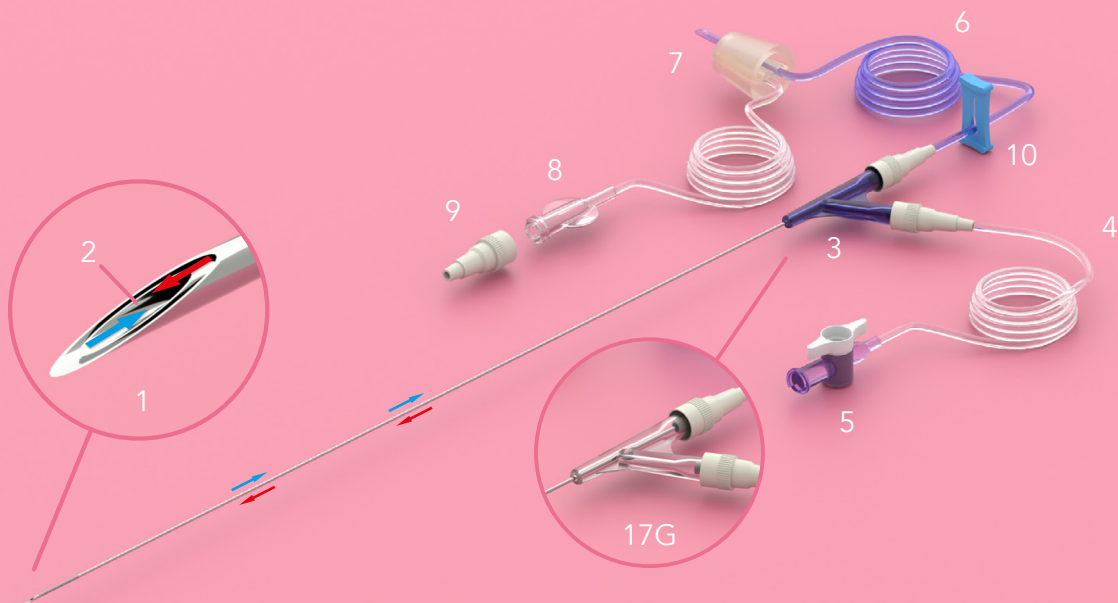

 Ago doppio lume a doppia via, la prima per aspirazione dell'ovocita (6) e la seconda per lavaggio follicolare (4), composto da:

- punta ecogena con affilatura atraumatica (1);
 - un lume interno dell'ago (2) per l'aspirazione ovocitaria prosegue attraverso un tubo azzurrato (6) che protegge l'ovocita dai raggi UV;
 - un tappo in gomma termoplastica medica (7) che si connette alla provetta;
 - un raccordo a Y (3), che funge anche da impugnatura, connette l'ago ai tubi di aspirazione e di lavaggio (provvisto di 2 raccordi universali);
 - un tubo in materiale trasparente lungo 100 cm (4) per il lavaggio follicolare termina con un rubinetto (5).
- Il lavaggio avviene attraverso il lume esterno dell'ago.
- una piastrina stringitubo (10) che regola l'aspirazione follicolare.

La connessione alla pompa di aspirazione è disponibile con prolunga di 80 cm che termina con un luer lock (8). Ulteriore adattatore universale (9) per una varia tipologia di connessioni.

 Double lumen, double-way needle, one way for oocyte aspiration (6), second way for follicle flushing (4) composed by:

- echogenic tip with atraumatic bevel (1);
 - oocyte aspiration is performed through the inner lumen of the needle (2) and continues through a blue-colored tubing set (6) protecting the oocyte from UV rays;
 - a thermoplastic medical rubber bung (7) connected to the test tube;
 - a Y connector (3), serving as finger grip, provided with two universal adapters connecting the needle to aspiration and washing tubing set;
 - follicle flushing is performed through a plastic tubing set 100 cm long (4) ending in a stopcock (5).
- Washing passes through the outer lumen of the needle.
- oocyte aspiration is controlled by a tubing clamp (10).

The connection to the suction pump is available through an extension tubing 80 cm long, ending in a luer-lock adapter (8). Additional universal adapter included (9) for a variety of connections.

Cod.	Colore - Color	Descrizione / Description	Pz/Scatola - Pcs/Box
770301/TL	Luer Lock Azzurro Light Blue Luer Lock	16 G x 330 mm	10 pz./pcs.
770308/TL	Luer Lock Trasparente Transparent Luer Lock	17 G x 330 mm	

