


Mininjector

Iniettori uterini per sonoisterografia
Sonohysterography uterine catheters

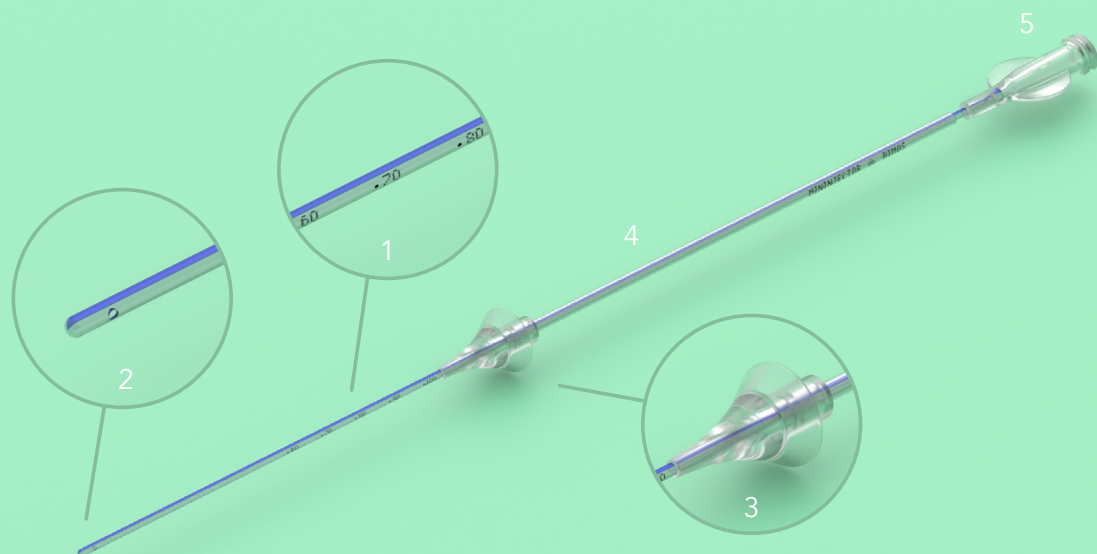

 Per l'iniezione di soluzione salina nelle procedure di sonoisterografia.

Il dispositivo è costituito da:

- un catetere trasparente centimetrato (1) per un posizionamento facilitato. La riga radiopaca blu lo rende ecograficamente visibile;
- una punta morbida con foro laterale (2);
- uno stopper conico (3) in morbido materiale plastico che ne assicura la tenuta sul collo dell'utero;
- un introduttore rigido (4) che facilita l'inserimento nel collo dell'utero;
- raccordo luer lock femmina (5) per la connessione della siringa.

 To infuse saline solution in sonohysterography diagnostic procedure.

The device is composed by:

- a clear catheter (1) with centimeter markings for easy positioning into the uterine cavity; a blue radiopaque stripe enhances visibility under ultrasound;
- a side-holed soft tip (2);
- a soft rubber stopper (3) for a guarantee of perfect sealing on the uterus neck;
- a rigid introducer (4) to advance into uterus neck;
- a female luer-lock adapter (5) for syringe connection.

Cod.	Descrizione / Description	Pz/Scatola - Pcs/Box
770314	Ø esterno 1,67 mm (5 Fr) O.D 1.67 mm (5 Fr)	15 pz./pcs.

Fornito in confezione singola - monouso - sterilizzazione a EtO - marcato CE 0051

Supplied in single packaging - disposable - EtO sterilization - CE 0051 marked

