


Monocyte-N

Ago per prelievo ovociti
Needle for oocyte aspiration

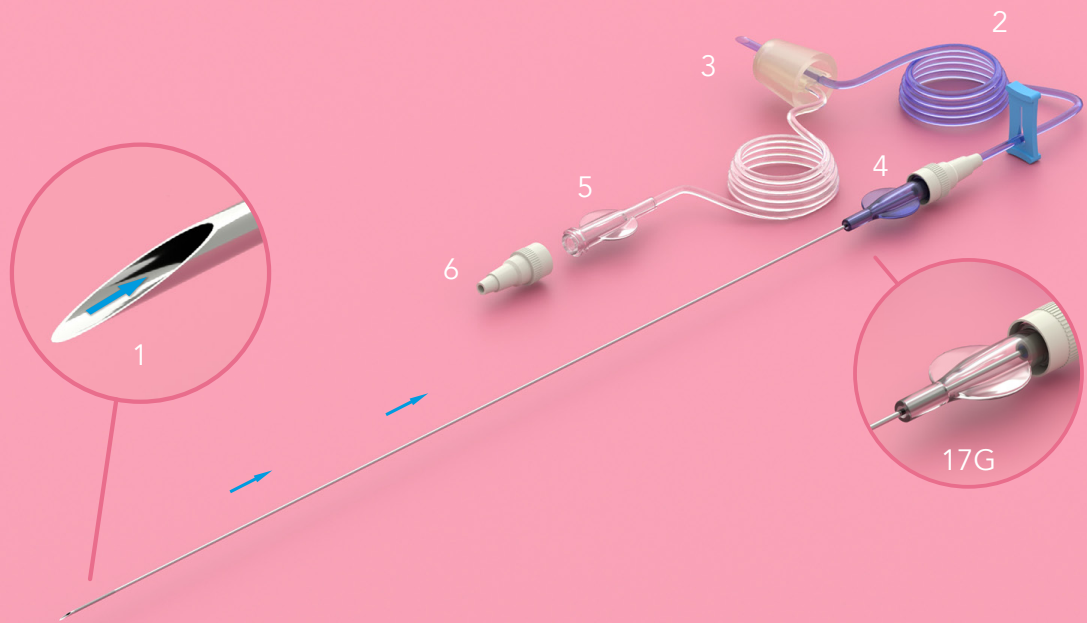

-  Ago monolume monovia composto da:
- punta ecogena con affilatura atraumatica (1);
 - tubo azzurrato (2) che protegge l'ovocita dai raggi UV durante l'aspirazione;
 - un tappo in gomma termoplastica medicale (3) che si connette alla provetta;
 - due alette laterali (4) in linea con l'affilatura della punta che permettono l'orientamento durante la fase di penetrazione.

La connessione alla pompa di aspirazione è disponibile con prolunga di 80 cm che termina con un raccordo luer lock (5). Ulteriore adattatore universale incluso (6) per una varia tipologia di connessioni.

-  Single lumen, one-way needle composed by:
- echogenic tip with atraumatic bevel (1);
 - blue-colored tubing set (2) that protects the oocyte from UV rays during aspiration;
 - thermoplastic medical rubber bung (3) connected to the test tube;
 - two orientation wings (4) aligned with tip bevel to orient the operator during penetration phase.

The connections to the suction pump is available through a 80 cm extension tubing ending in a luer lock adapter (5). Additional universal adapter included (6) for a variety of connections.

Cod.	Colore - Color	Descrizione / Description	Pz/Scatola - Pcs/Box
770300/TL	Luer Lock Azzurro Light Blue Luer Lock	16 G x 330 mm	10 pz./pcs.
770307/TL	Luer Lock Trasparente Transparent Luer Lock	17 G x 330 mm	

Fornito in confezione singola - monouso - sterilizzazione a EtO - marcato CE 0051

Supplied in single packaging - disposable - EtO sterilization - CE 0051 marked

